

# ergotron®

*moving you forward*


## **StyleView 42 Medical Cart**

- Easy to maneuver: Moves easily across carpet or hard surface floors with a small footprint that fits in space-constrained areas
- Enhances well-being and productivity: Caregivers can choose their preferred height for ergonomic typing and charting as they sit or stand
- Promotes patient and caregiver safety: Entire power system and cart are certified to UL/EN/IEC 60601-1 and comply with EMI limits of FCC Part 15 Class A
- Long battery life: LiFe battery option helps improve uptime with an ultrafast recharge and long runtime that can last up to 12 hours, depending on the number of connected devices

- Quick integration: Ships fully assembled

### *Electric lift height adjustment option*


- Accessibility for all: Reduces user forces with heavy cart loads or for caregivers with physical limitations
- Adjustment at the push of a button: Fits a broad range of users 5'0" to 6'1" (152 to 185 cm) with height adjustment that moves the worksurface and keyboard

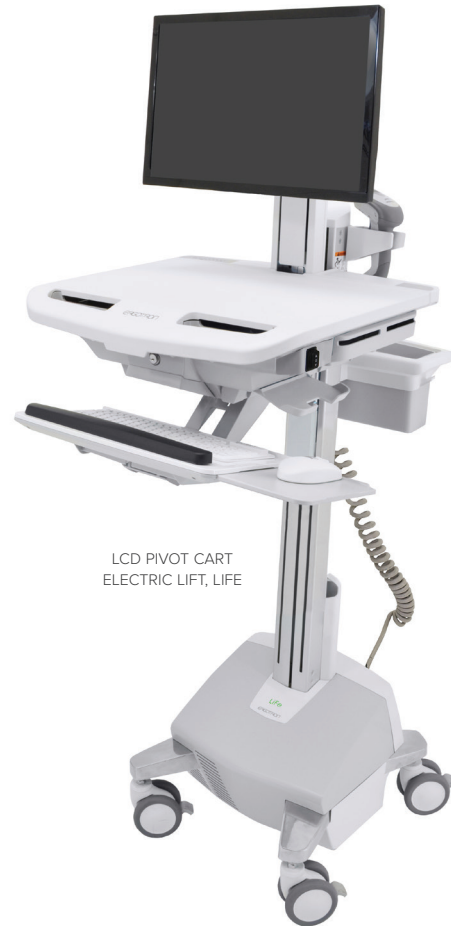

LAPTOP



LCD ARM

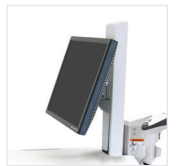

LCD PIVOT

## VERSATILE POWERED CART SUPPORTS MULTIPLE WORKFLOWS

## STYLEVIEW™

StyleView Medical Carts bring comfort and flexibility to any workflow by adapting to meet the needs of every caregiver. Built on decades of experience, these compact carts offer top-notch ergonomics and versatility. The open architecture design, broad portfolio of accessories and quality construction makes StyleView a long-term investment that evolves with your organization.

Customize your fleet to fit your workflows with:

*Three display mount options*

*Lithium Iron Phosphate (LiFe) or Sealed Lead Acid (SLA) battery power*

*Proprietary Constant Force™ (CF) manual height adjustment or electric lift (eLift)*

# STYLEVIEW™

## PRODUCTS INCLUDE:

LCD Pivot or LCD Arm mounting options, antimicrobial worksurface, CPU lock, secure CPU compartment, four-port USB hub, scanner holder bracket, mouse holder, wrist rest, negative-tilt keyboard tray, left/right mouse tray, rear handle, storage basket, cable storage, four casters, eKinnex™ fleet management software.

- 1 Convenient scanner holder bracket supports smooth workflows
- 2 Large worksurface and wrist rest are treated with antimicrobial to support infection control
- 3 Negative-tilt keyboard tray, wrist rest and mouse holder support comfortable documentation; four amber LEDs provide visibility of the keyboard tray in low-light care areas
- 4 Rear storage basket keeps task-related items close for easy access
- 5 T-slot channel interface on the column allows you to integrate accessories
- 6 Four 4" (10.2 cm) lockable casters, including one with a tracking function to easily steer while moving the cart
- 7 Compact base fits into crowded healthcare spaces



LCD PIVOT CART  
CONSTANT FORCE, LiFe

Durable, easy-to-clean exterior composed of aluminum, high-grade plastic and zinc-plated steel



T-SLOT WIPES HOLDER  
98-450



4 DRAWERS (1 PRIMARY TRIPLE AND  
1 SUPPLEMENTAL SINGLE)



9 DRAWERS

Additional accessories available at [healthcare.ergotron.com](http://healthcare.ergotron.com). Drawers can be added or removed based on changing needs. Please note, adding drawers may impact the total height adjustment.

	CF Laptop, LiFe	CF LCD Arm, LiFe	CF LCD Pivot, LiFe	eLift LCD Arm, LiFe	eLift LCD Pivot, LiFe
Part # (white)	SV42-6102	SV42-6202	SV42-6302	Americas, APAC: SV42-7202	Americas, APAC: SV42-7302
Typical LCD size	≤ 17.3"	≤ 24"			
LCD SIZE IS APPROXIMATION. DIMENSION IS DIAGONAL MEASUREMENT OF SCREEN. LCD SIZE MAY BE EXCEEDED AS LONG AS SCREEN WEIGHT DOES NOT EXCEED MAXIMUM WEIGHT CAPACITY OF THE MOUNT					
Weight Capacity*	≤ 18 lbs (8.2 kg)	≤ 31 lbs (14.1 kg)			
Lift*	20" (51 cm)	27" (68.6 cm)	25" (63.5 cm)	27" (68.6 cm)	25" (63.5 cm)
Tilt*	—	65°/+60°/-5°	25°/+20°/-5°	65°/+60°/-5°	25°/+20°/-5°
Pan*	—	180°	24°	180°	24°
Rotation*	—	360°	90°	360°	90°
VESA	—	MIS-D (100 x 100 mm & 75 x 75 mm)			
Cart Weight	93 lbs (42.3 kg)	97 lbs (44.1 kg)			
Shipping Dimensions	56"H x 34"W x 24.5"D (142 x 86 x 62 cm)				
Shipping Weight	127 lbs (57.7 kg)	131 lbs (59.5 kg)			
Warranty	5-year warranty for mechanical components; 3-year warranty on LiFe batteries, electric lift and switch; 2-year warranty on power modules and power coiled cords. See Ergotron product warranty terms & conditions for more detail.*				

\*See website

SLA power options available in some regions, contact Ergotron for more details.

Download additional resources at [ergotron.com](http://ergotron.com).

For more information:

USA & Canada: 800.888.8458 / +1.651.681.7600 / [insidesales@ergotron.com](mailto:insidesales@ergotron.com)

EMEA: +31.33.45.45.600 / [info.eu@ergotron.com](mailto:info.eu@ergotron.com)

APAC: [apaccustomerservice@ergotron.com](mailto:apaccustomerservice@ergotron.com)

LATAM: +1.800.681.1820 / [info.latam@ergotron.com](mailto:info.latam@ergotron.com)

Custom: [custom@ergotron.com](mailto:custom@ergotron.com)

© Ergotron, Inc. 11.09.2023

Literature made in the U.S.A.

Content subject to change.

Ergotron devices are not intended to cure, treat, mitigate or prevent any disease.

Patent information available at [www.ergotron.com/patents](http://www.ergotron.com/patents).

